


Vertical Autoclave

This vertical autoclave is especially designed for liquid and food sterilization in laboratory, pharmaceutical and food factory.

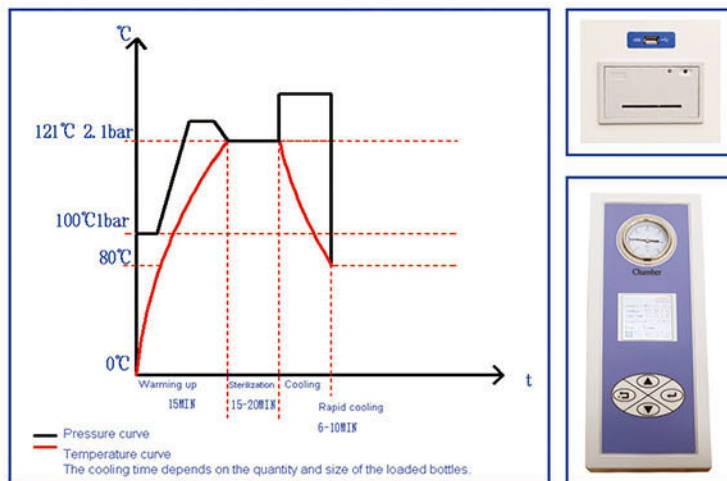

Features:

- * Standard build-in printer.
- * Deionized water interface for supplying steam generator.
- * F0 value display and printed.
- * Rapid water cooling system for efficient and safe sterilization of liquids.
- * With a quick maintenance window, the electrical parts can be repaired without removing the cover.
- * Microcomputer control. Fully automatic control water injection, temperature rise, sterilization, exhaust and dry.
- * LCD screen display shows pressure, temperature, time, running status, error codes and operation curve etc.
- * Mobile probe detects the internal temperature of liquid directly, to ensure the sterilization effect of liquid program.

Technical Parameters:

Model	BKQ-B50L	BKQ-B75L
Capacity	50L	75L
Chamber Size(mm)	φ386*490	φ386*670
Chamber Material	SUS304	
Designed Pressure	0.28Mpa	
Designed Temp	150°C	
Working Pressure	0.22MPa	
Working Temp.	80°C~136°C	
Temp. Display Accuracy	0.1°C	
Noise	≤65dB	
Consumption	5.5KW	
Power Supply	AC110/220V±10%, 50/60Hz	
SS Basket	1 pc	
External Size(W*D*H)mm	700*610*1100	
Package Size(W*D*H)mm	800*715*1270	
Gross Weight(kg)	140	147