



Configuration
UV Light;
Fluorescent Light;
Adjustable Spotlight;
Stainless Steel Hanger;
"UV, FAN, LIGHT,
PULVERIZER" Control.




Water Tap & Sink
Stainless steel;
Whole-structure.



Suction Inlet
Suction inlet adopts louver-type.
The upper mouthing structure
ensures above blower suction,
avoiding air back-streaming.



**Flex Magnifier
Light(Optional)**

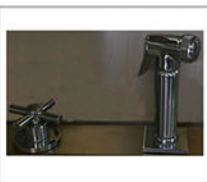

Auto Flushing Design
Automatic flushing device, slightly
right-side slope worktable, easy to
waste discharge and cleaning.



Pulverizer (Optional)
Completely smash the waste
during test, to avoid blocking
drainage pipes.

APPLICATION

Pathological sampling bench is widely used in hospital pathology department, pathology laboratory, etc. Reasonable ventilation system protects operator away from the harmful gas produced by the formalin during sampling. The hot & cold water system ensures it can adapt work in different climates.

ADVANTAGE

- 1. Odor-remove function.
- 2. Single water sink, two-person operations.
- 3. Exhaust air system: High-power suction device, low noise.
- 4. Optional water heater is available, realizing hot & cold water adjustment.
- 5. The base cabinet is inward-concave type, convenient for sitting operation.
- 6. Optional pulverizer can completely smash the waste during test, to avoid blocking drainage pipes.
- 7. Automatic flushing device, slightly right-side inclined worktable, easy for waste discharge and cleaning.

TECHNICAL PARAMETERS

Model		QCT-1000	QCT-1500	QCT-1800
External Size(W*D*H)		1000*750*1950mm	1500*750*1950mm	1800*750*1950mm
Working Area Size(W*D*H)		895*430*900mm	1395*430*900mm	1695*430*900mm
Work Surface Height		900mm		
Consumption		250W	300W	300W
Optional Water Heater		1500W		
Power Supply		AC220V±10%, 50/60Hz; 110V±10%, 60Hz		
Exhaust Airflow Volume		500~600m³/h		
Material	Main Body	304 stainless steel		
	Work Table	304 stainless steel		
Standard Accessories		Spotlight; UV lamp; Fluorescent lamp; Exhaust system; 4m exhaust duct; Pipe strap		
		Auto flushing device; Stainless steel hanger; Adjustable hot and cold water tap; 180mm deep stainless steel water sink, funnel with filter. Active carbon filter		
Optional Accessories		Pulverizer, Electrical water heater, Automatic induction hand dryer, Eye washer, Camera, Flex magnifier light, laptop bracket		
Package Size(W*D*H)		1160*890*2170mm	1660*890*2170mm	1960*890*2170mm
Gross Weight		250kg	300kg	340kg